
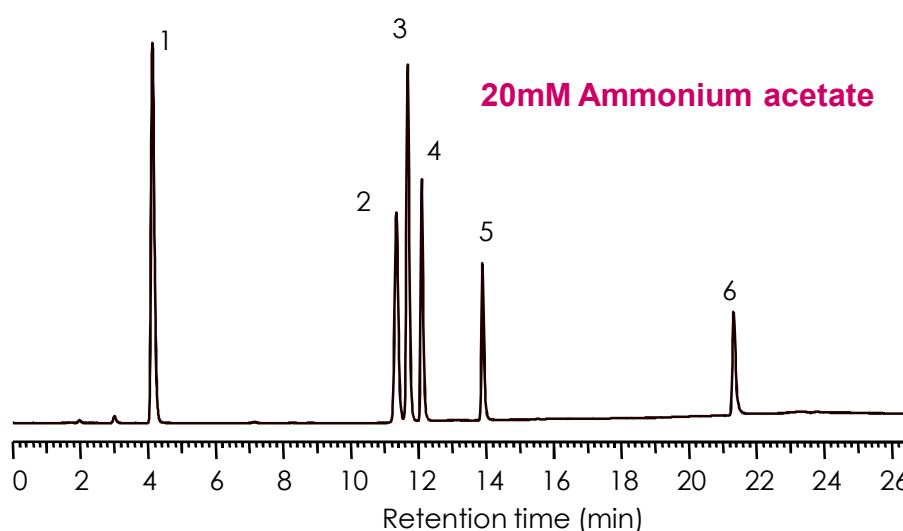
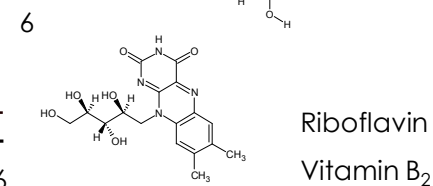
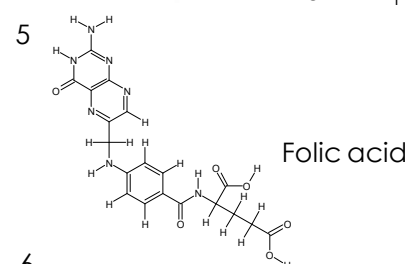
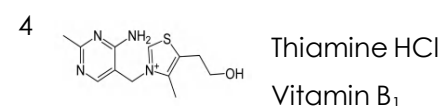
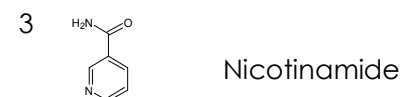
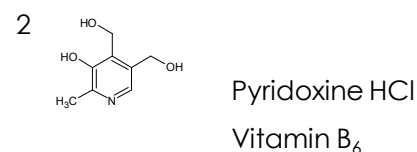


水溶性ビタミン類

Water-soluble vitamins

Sunniest RP-AQUA 5 μ m, 150 x 4.6 mm i.d.

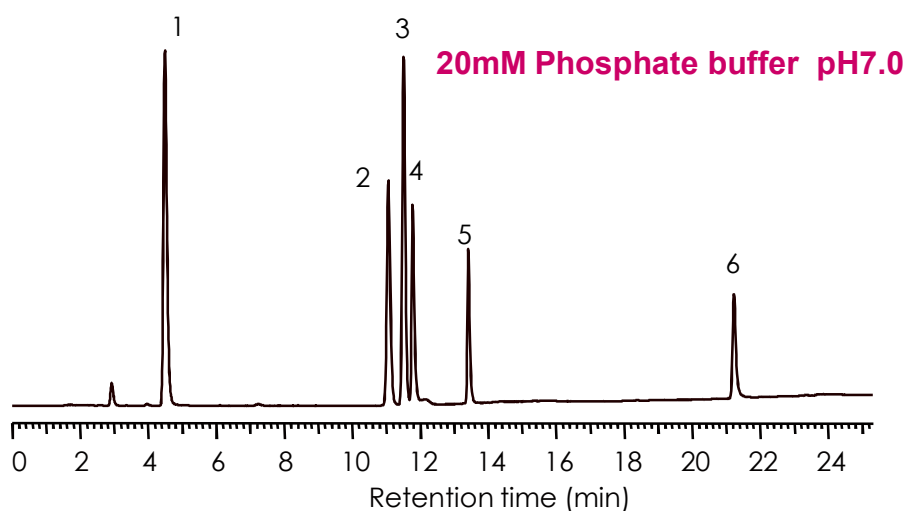

Column: Sunniest RP-AQUA 5 μ m, 150 x 4.6 mm i.d.

Mobile phase: A) 20mM Phosphate buffer pH7.0 or 20mM Ammonium acetate

B) Acetonitrile/ A solution (20:80)

Time (min)	0	5	20	25
%B	0	0	100	100

Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@250 nm

Injection volume: 10 μ L

★LC/MS用の移動相（酢酸アンモニウム）でもリン酸緩衝液と同じ分離