


BSAの消化物の分離

SunShell C18-16 2.6 μ m, 100 x 2.1 mm i.d.

BSA Digest using SunShell C18-16

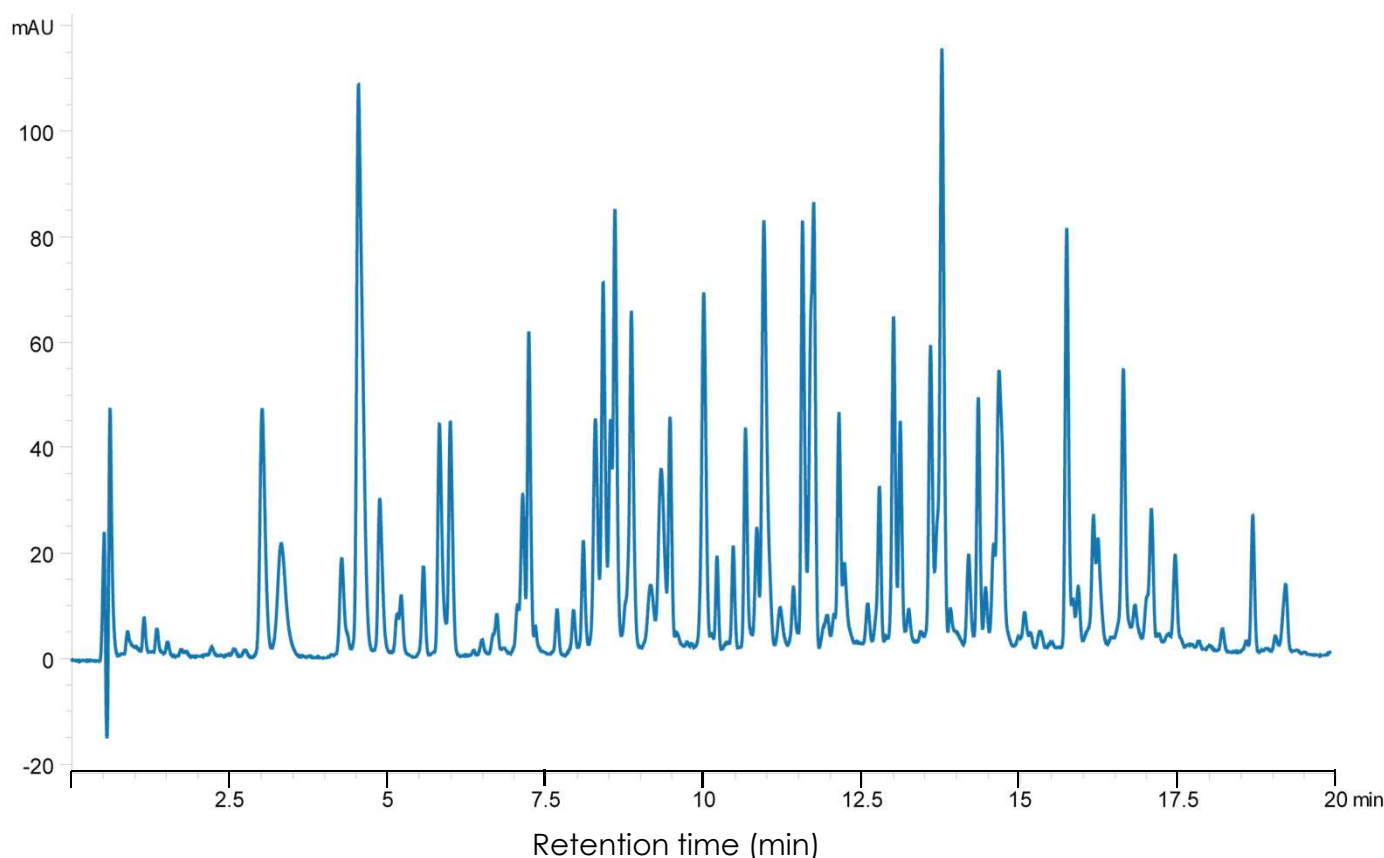

Column: SunShell C18-16, 2.6 μ m (16 nm) 100 x 2.1 mm,

Mobile phase: A) 0.1% TFA in water

B) 0.1% TFA in acetonitrile

Gradient program:

Time(min)	0	1	21
%B	10	10	40

Flow rate: 0.50 mL/min ,

Temperature: 40 °C

Detection: UV@214 nm,

Injection volume: 10 μ L,

Concentration: 0.01 μ g/ μ L

Sample:

1. Bovine Serum Albumin Digest