
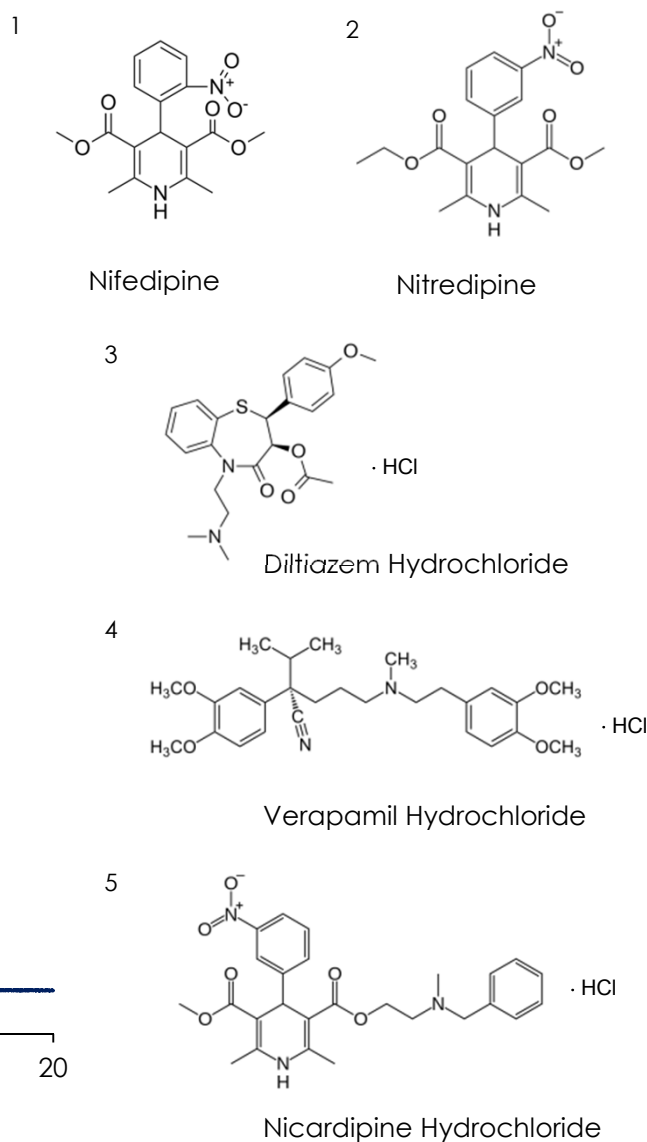


Ca拮抗剤の分離 (Sunniest C18)

Calcium antagonists

Sunniest C18 5 μ m, 150 x 4.6 mm i.d.

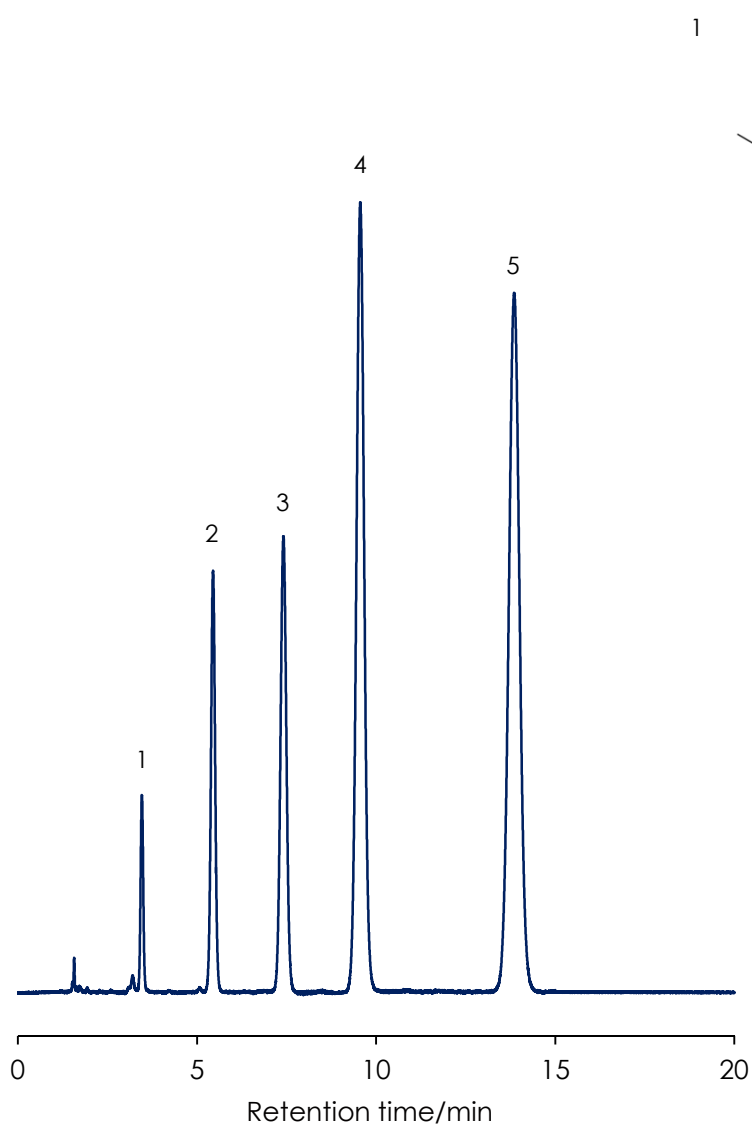

Column: Sunniest C18, 5 μ m 150 x 4.6 mm

Mobile phase: Methanol:25 mM phosphate buffer pH7=70:30

Flow rate: 1.0 mL / min

Temperature: 40 °C

Detection: UV@230 nm

Sample: 1 = Nifedipine, 2 = Nitredipine, 3 = Diltiazem Hydrochloride, 4 = Verapamil Hydrochloride,

5 = Nicardipine Hydrochloride