


## 水溶性ビタミン類 (3)

### Water-soluble vitamins (3)

SunShell RP-AQUA 2.6  $\mu$ m, 150 x 4.6 mm i.d.

#### SunShell RP-AQUA

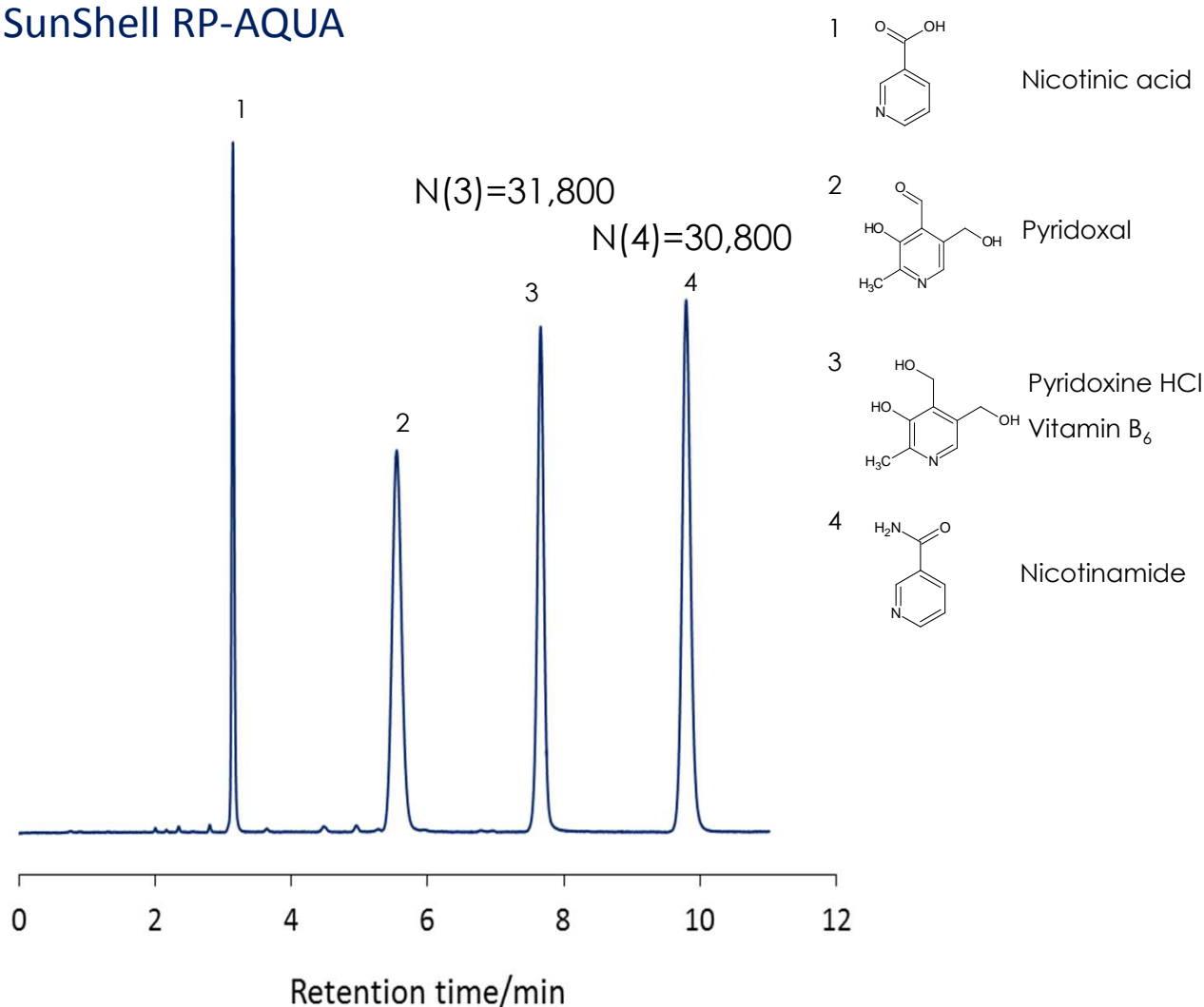

Column: SunShell RP-AQUA 2.6  $\mu$ m, 150 x 4.6 mm i.d.

Mobile phase: 40mM Phosphate buffer pH6.8

Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@250 nm

Sample: 1 = Nicotinic acid, 2 = Pyridoxal, 3 = Pyridoxine HCl, 4 = Nicotinamide

UHPLC: X-LC (Jasco)