
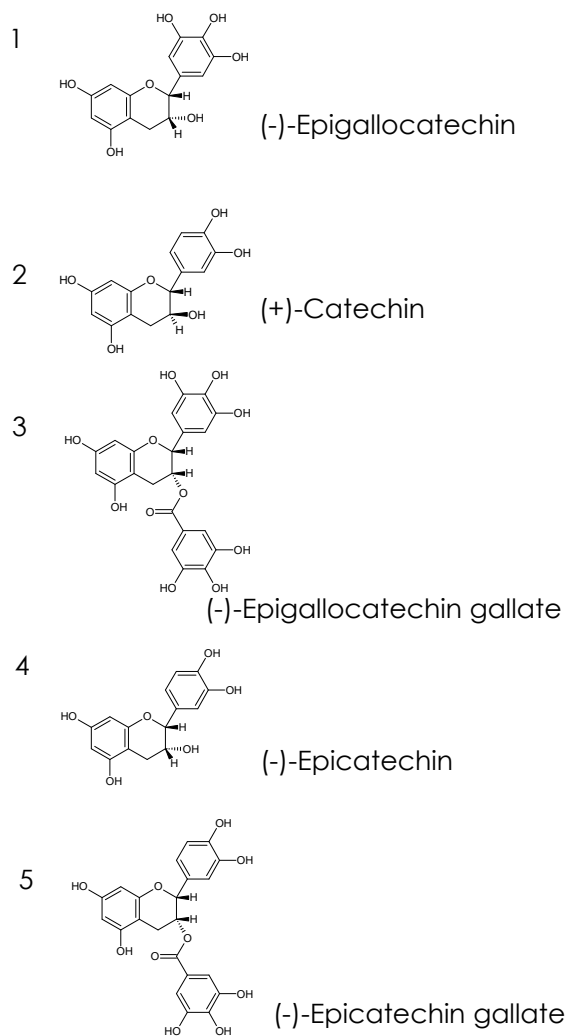


カテキン類 (2)

Catechins (2)

Sunniest RP-AQUA 5 μm, 150 x 4.6 mm i.d.

Sample: Oolong tea

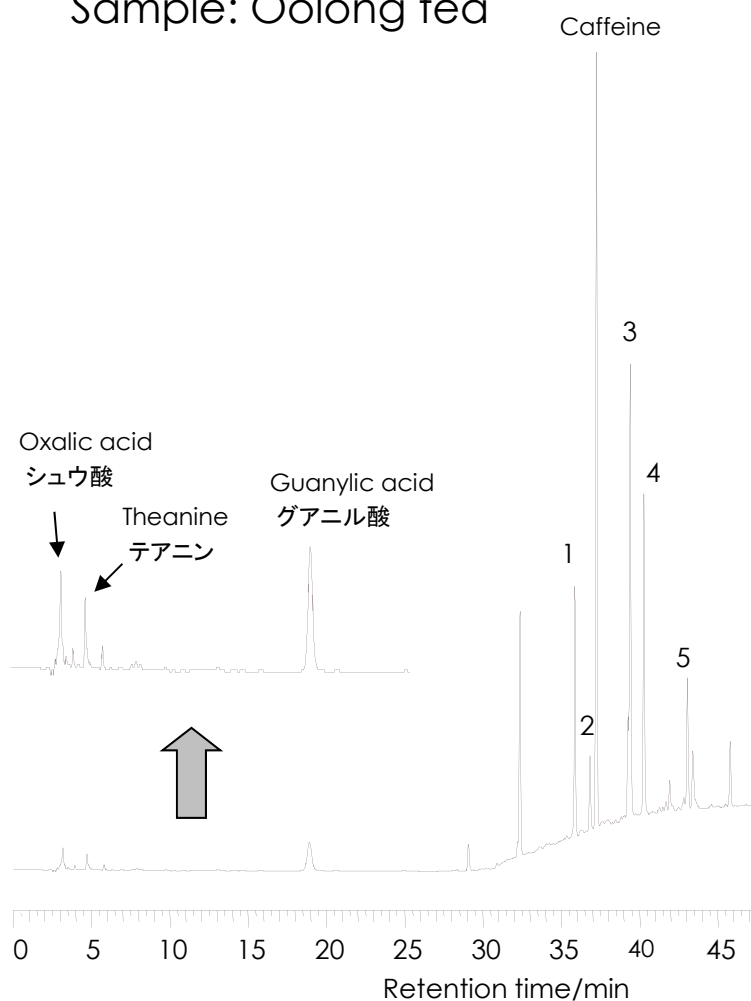

Column: Sunniest RP-AQUA 5 μm, 250 x 4.6 mm i.d.

Mobile phase: A) H₂O/H₃PO₄ (99.9:0.1)

B) CH₃CN

Time (min)	0	20	60
%B	0	0	50

Flow rate: 1.0 mL/min

Temperature: 40 °C

Detection: UV@210 nm

Injection volume: 10 μL