



多糖(ガラクトン)の分離

SunBridge C18 5 μ m, 150 x 4.6 mm i.d.

Polysaccharides (galactan)

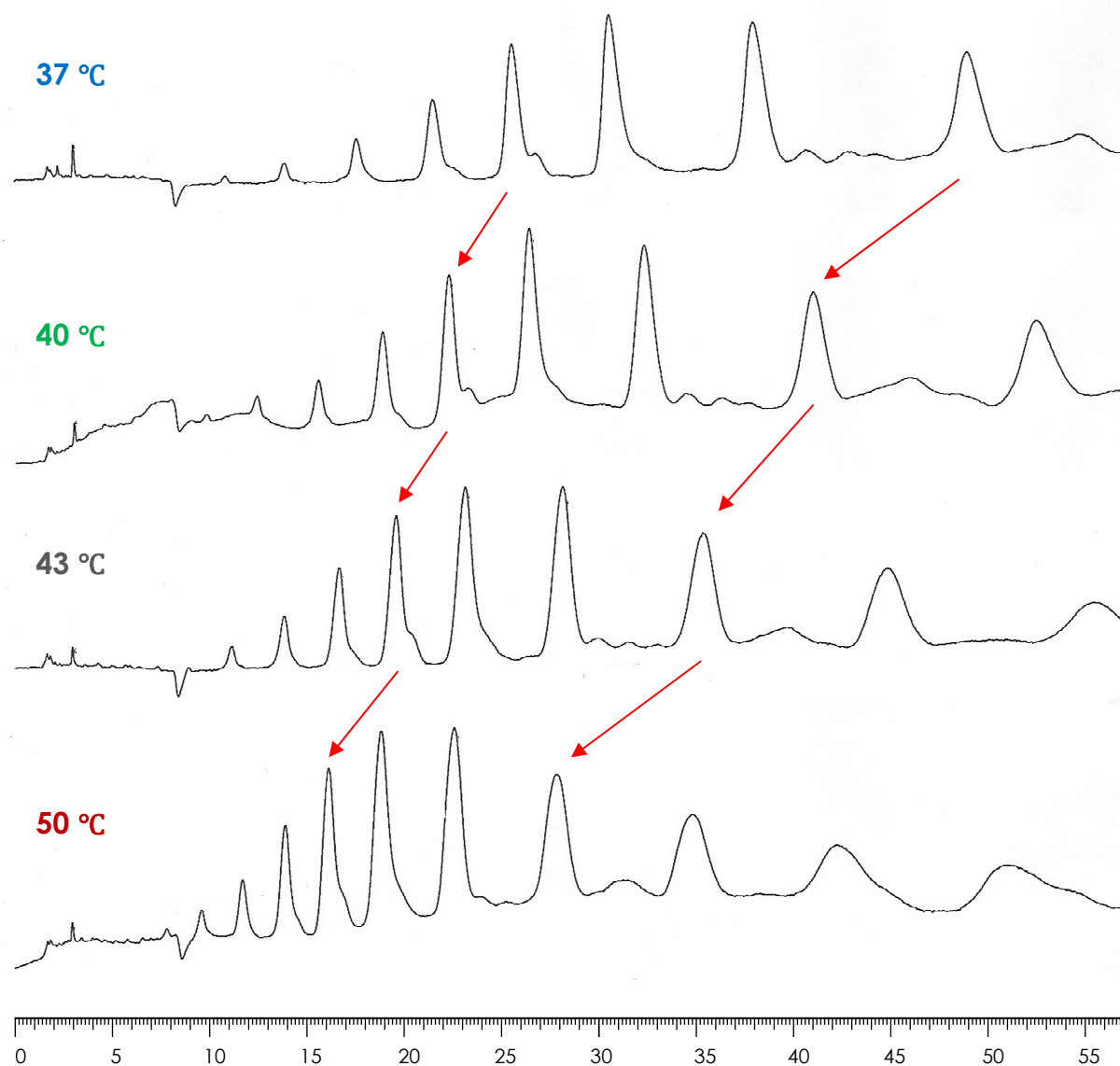

Column: SunBridge C18 5 μ m 150 x 4.6 mm
Mobile phase: Water
Flow rate: 1.0 mL / min
Temperature: 37 °C, 40 °C, 43 °C, 50 °C
Injection volume: 50 μ L
Detection: RID (Shodex RI-101)
Samples: Sample: Galactan oligomers (DP $\approx$ 16–25)

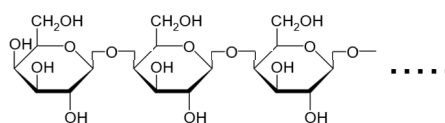

Galactan: composed of galactose units

提供: 富山県立大学大学院工学研究科 酵素化学工学講座 加藤 康夫先生